BEAUTIFLY

Lumi
Blush
series

Depilator laserowy
Laser hair remover

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beautifly.eu

@beautifly.eu

beautifly.eu

/beautiflyeu


@beautifly.eu

f1

f2

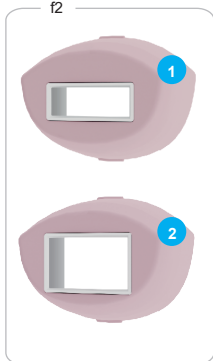f3

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- www.beautifly.eu
- www.youtube.com/beautiflyeu

Input power:	36W
Rated voltage:	100-240V
Rated frequency range:	50Hz - 60Hz
Storage temperature:	from 0oC to 40oC
Working temperature:	from 5oC to 35oC
Humidity during operation:	from 30% to 85%

Security:

Before using the product, read this manual and follow the recommendations contained therein when using the device.

General Description [page 2- f1]:

1. Epilator window
2. Screen
3. On/Off switch
4. Skin assessment/test button
5. Pulse button
6. Ventilation
7. Power input

Function [page 2- f2]:

B-Lumi has two functions and separate lamps for each of them.

1. IPL hair removal
2. IPL skin revitalization

Auto Pulse Mode:

Press the pulse button for 2 seconds, the screen will show (double lightning icon), at the same time the device will switch to automatic mode. The device will flash automatically when the pulse window is on your skin. When the device is in automatic mode, press the pulse button for 2 seconds, the device will switch to manual mode shown (single lightning icon).

Contraindications:

1. The device is not intended for children under 15 years of age. Young people aged 15–18 may use the device with the consent and/or supervision of their parents or guardians.
2. Do not use if the skin in the area to be treated is tanned or red. Using the device on tanned or red skin may cause temporary or even permanent discoloration of the skin.
3. Do not use if you are pregnant, trying to get pregnant, think you may be pregnant or are breastfeeding. The B-Lumi device has not been tested for this purpose and therefore should not be used in these cases. Hormonal changes that occur during pregnancy or breastfeeding may increase the risk of skin damage or skin sensitivity.
4. Do not use the device if you suffer from diseases such as epilepsy, epileptic seizures or migraine. Energy flashes can cause severe headaches and migraines.
5. Do not use on areas of skin covered with tattoos, dark brown or black spots (such as dark freckles, moles, or warts). Tattoos or dark spots can absorb too much light energy, which can cause skin swelling, temporary or permanent blisters, discoloration, or skin whitening.
6. Do not use the device if you have a tendency to develop keloid scars.
7. Do not use on the scalp, ears, neck, nipples or mucous membranes (labia minora, vaginal and anal areas). In those parts of the bikini area where the skin is more pigmented, a lower light intensity level should be used. Use in these areas may cause skin swelling, blistering, or temporary discoloration or whitening of the skin.
8. Do not use if the treatment area has any active skin inflammation, including acne, psoriasis, eczema, burns, herpes, abrasions or rash. These factors can make the skin more sensitive to light and make the treatment uncomfortable or cause swelling, blisters or infections on the skin. Once your skin has healed sufficiently, consult your doctor before using the device.
9. Do not use the device if you are taking medications that make your skin sensitive to light, including non-steroidal anti-inflammatory agents (e.g. aspirin, ibuprofen, cetaminophen), tetracyclines, phenothiazines, thiazide, diuretics, sulfonyleureas, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha Hydroxy Acids (AHAS), Beta-Hydroxy Acids (BHAS), Retin-A, Accutane and/or topical retinoids.
10. Do not use if you have a disease associated with photosensitivity, e.g. porphyria, polymorphic photoeruption, solar urticaria, lupus, etc. These conditions make the skin more sensitive to light and the procedure may be unpleasant or causing skin swelling, blisters or infections.

Action Procedure:

To get the best results, we suggest that you read the user manual before first use including all safety instructions and warnings for correct and safe use of the device. The Lumi hair removal device has integrated filters that protect the skin against harmful radiation and burns during the hair removal treatment. The rate of hair growth varies, depending on age, metabolism and other factors. Growth hair is usually divided into three stages.

Growing phase:

Hair actively grows along with the bulbs. At this stage, melanin concentration is the highest. Melanin determines the formation of hair color. The concentration of melanin is also a decisive factor intensive pulsed light technology activities. Only hair in the growing phase will react to it light hair removal process.

Degeneration phase**Stopping hair growth and rooting the hair bulb.****Reconstruction phase**

Mature hair detaches from the bulbs and falls out. The bulb is dormant until the next cycle hair growth.

In the first stage of hair removal with the Lumi device, light energy is used to operate directly onto the hair roots below the skin. Melanin in the skin and hair follicles absorb impulse energy of light emitted by the device. The darker the hair color, the more light energy is absorbed, and this process stimulates the bulbs to enter the dormant phase.

In the second stage, during the depilation period, the hair will fall out naturally and the regeneration will continue inhibited, which may take 1-2 weeks. Only hair in the growing phase will be removed. To effectively remove unwanted hair, repeat the depilation process once a week for a period of time two months.

Handling Warning:

1. The product will not turn on without inserting the lamp, please insert the lamp before connecting to the power. When the device is operating, do not cover the air intake and ventilation openings to prevent this had a negative impact on heat release.
2. Wear safety glasses when operating the device.
3. Clean the area of skin being treated before use. It is not allowed to use any cosmetics on the skin immediately before the treatment. This may cause destruction device and influence the effect of the treatment.
4. Adjust your energy level from low to high. The higher the energy level, the better result. Immediately after the treatment, skin redness may occur.
5. After phototherapy treatment, you should use a moisturizing preparation and preparations with high protection sunscreen.
6. After finishing the treatment, use a damp cotton swab or cloth to clean the lamp from hair residues, sebum or sebum.
7. It is recommended that the device and accessories be stored away from water and moisture. Store the device in a ventilated and dry environment and keep it away from kids.
8. It is prohibited to place metal objects in the heat emission hole. Otherwise there is a risk of electric shock.

Tip Installation:

Use your thumb and index finger to hold the strap-like bulge on the sides up and down down, push into the slot to complete installation.

Removing the Tip:

Use your thumb and index finger to hold the strap-like bulge on the sides of the up and down, pull along the slot.

Take the Skin Sensitivity Test:

If you are using the device for the first time or on a new area of the body, select a small area in near the area you intend to treat to perform a skin sensitivity test. Use level of the lowest device. Follow the instructions below to perform the procedure, but emit no more than 2-3 pulses. After performing the skin test, wait 60 minutes and check if the skin does not react. If after 60 minutes the skin appears normal and there is no redness or other skin reaction, you are ready for the full treatment.

Choose the Right Intensity:

Choose the appropriate intensity depending on your skin tolerance and increase it gradually from one up to five degrees. If the intensity is too weak and there are no skin reactions, you can gradually increase the intensity level. Do not flash pulses in one place on the skin multiple times. You this creates too much heat on the skin, which could cause the skin to swell, blister or become temporary or permanent changes in the color of your skin. After each flash, slowly move the device to adjacent place on the skin. Wait at least before returning to the previous area of skin at least 20 seconds.

Hair removal:

Application: Mustache, armpit hair, body and leg hair, hairline hair on the forehead and bikini area. The effect varies depending on skin color. The device is not suitable for whites and bright.

Usually it is about 7-9 treatments. Perform the first 3 treatments once a week. For the next 4-9 treatments, perform one therapy every 2-3 weeks. To maintain the effect of hair removal and to prevent the hair follicles from reactivating, repeat the treatments regularly to ensure that the skin remains smooth. If treatments are not continued, the hair follicles will become active again and the hair will begin to grow back. To maintain the results, simply repeat the treatments every 4-8 weeks to ensure that the hair follicles do not become activated and hair regrowth remains inhibited. The intervals between treatments may vary depending on the individual rate of hair regrowth.

Table comparing skin and hair color

How often to perform hair removal treatment:

✓ applicable

✗ not applicable

Skin color \ Hair color	White	Ivory tint	Beige	Light	Brown	Dark brown and black
Black	✓	✓	✓	✓	✓	✗
Dark brown	✓	✓	✓	✓	✓	✗
Brown	✓	✓	✓	✓	✓	✗
Golden	✓	✓	✓	✓	✓	✗
Red	✗	✗	✗	✗	✗	✗
White	✗	✗	✗	✗	✗	✗
Gray	✗	✗	✗	✗	✗	✗

How to Perform a Hair Removal Treatment:

1. Before using the Lumi device, prepare your skin by shaving the selected area to be treated. The skin should also be washed.
2. Wear safety glasses to avoid IPL light stimulating your eyes.
3. Select the appropriate intensity (see "select intensity") and start at a low level.
4. Let the device window touch the selected areas vertically and treat the areas one by one. Don't repeat the same place over and over or skip any area.
5. Start on the initial area and repeat step 4, 2-3 times.
6. After the treatment, it is recommended to use a moisturizing preparation after depilation.
7. Turn off the device and use a damp cloth or cotton swab to clean the treatment window.
8. Place the device and accessories in a cosmetic bag or box and store them in a safe place.

Skin photorejuvenation:

For people who want to visibly rejuvenate the skin and refresh its appearance, evening out the color and reducing troublesome discolorations, we offer photorejuvenation treatments with the modern Lumi cosmetic laser. Lumi laser therapy allows for significant smoothing of the skin and complete removal of age spots and erythema.

The aim of rejuvenating selected areas is to use short and safe light pulses directed towards the skin surface. The light wave stimulates skin cells to produce collagen and elastin, which accelerate skin regeneration. The therapeutic effects of IPL radiation include improved blood supply and oxygenation of the skin, an increase in the number of collagen and elastic fibers, better moisture absorption, and regulation of the pigmentation and keratinization process.

Effects:

- Elimination of symptoms of skin aging
- Reconstruction of collagen fibers
- Firm and elastic skin
- Reduction of pores
- Smoothing out wrinkles
- Improving skin hydration
- Lightening discolorations

How Often to Perform a Photorejuvenation Treatment:

To obtain the desired effects, a series of treatments must be performed. It is recommended to carry out a series 8 treatments, one treatment each week.

How to Perform a Photorejuvenation Treatment:

1. Before starting the treatment, you should thoroughly cleanse your skin of makeup and any cosmetics. You should also wipe the lamp.
2. Wear safety glasses to avoid IPL light stimulating your eyes.
3. Select the appropriate intensity (see "select intensity") and start at a low level.
4. Let the device window touch the selected areas vertically and treat the areas one by one. Don't repeat the same place over and over or skip any area.
5. Start on the initial area and repeat step 4, 2-3 times. After the treatment, it is recommended to use a strongly moisturizing preparation.
6. Turn off the device and use a damp cloth or cotton swab to clean the treatment window.
7. Place your device and accessories in a case or box and store them in a safe place.

More information about the company and products [page 2 - f3]:

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled.

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/UE of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/UE of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/UE of the European Parliament and of the Council as regards the list of restricted substances.

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