Drogi Kliencie!

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
HX701/HX702 | REV.A | 03.2023



BATTERY REMOVAL


1. The process is irreversible.
2. Please follow necessary safety procedures and protect your eyes, hands, and fingers.

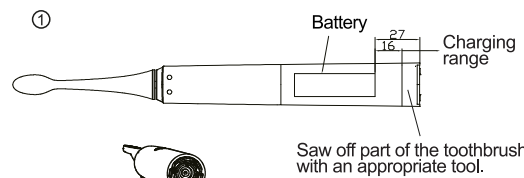Saw off part of the toothbrush with an appropriate tool.

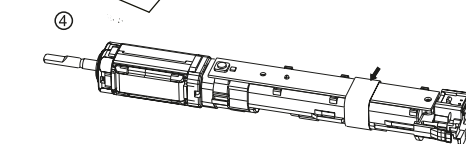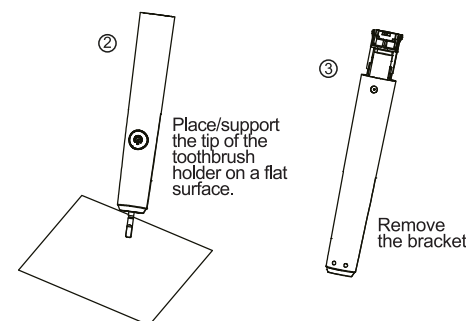After removing the bracket, a transparent tape will be visible, and the place marked with an arrow (the transparent tape) should be cut open.

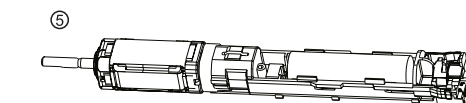After cutting the transparent tape, the battery can be removed from the support bracket.

UTILIZATION

INFORMATION FOR USERS ON HOW TO HANDLE WASTE ELECTRICAL AND ELECTRONIC EQUIPMENT

The crossed-out wheeled waste container symbol on products indicates that waste electrical and electronic equipment may not be disposed of with municipal waste. Disposing of electrical waste in a container poses a threat to the environment. By segregating and submitting waste electrical and electronic equipment for treatment, recovery, recycling and disposal, you protect the environment from pollution and contamination, help reduce the use of natural resources and reduce the production costs of new equipment. For details on the nearest collection point, please contact your dealer or company:

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USING THE TOOTHBRUSH

1. Connect the charging port to the power outlet and start charging. When the brush is fully charged, disconnect it from the cable before use.
2. Insert the brush tip and make sure it is firmly attached. The gap between the tip and the brush, approximately 0.5 mm wide, is due to the design of the device. **DO NOT FORCE THE TIP IN.**
3. Soak the bristle with water and apply an appropriate amount of toothpaste.
4. Bring the toothbrush close to your teeth.
5. Switch on the toothbrush and gently swipe its tip on the teeth, devoting 10 seconds to each area of the oral cavity, which should cover the front, rear and top tooth areas.
6. When using the mode switch function for the first time, push the start button after activating the device to switch between available modes. Continuous pressing of this button for 3 seconds will automatically turn off the toothbrush.
7. Press the multifunctional button to select from among five modes depending on individual preferences: whitening, cleaning, sensitive, polishing and massage.
8. When the brushing time exceeds 30 seconds for a single oral area, the brush will automatically remind you of the need to change its position by means of a short pause. Two minutes after it is started, the brush will be automatically switched off.

CHARGING INSTRUCTIONS


1. When the battery is low, the LED indicator will blink 10 times. NOTE: The charging indicator will not light up when the battery level is lower than 10%.
2. Once the device is connected to the power source, the battery charging process will start and the LED indicator will light red. When the battery is fully charged, the indicator will change to green.
3. When the toothbrush is used for the first time or when it has not been used for a long time, the indicator may not light up for the first 5 minutes - this is completely normal and due to its low battery level.

PRODUCT SPECIFICATIONS


Model	HX701/HX702
Numbers of programs	5
Input voltage	5 V
Input power	1.5 W
Battery type	lithium battery 700mAh
Product material	ABS/ TPE
Blisters material	Dupont Nylon
Size	26.9 x 26.4 x 247 mm
Weight	102 g
Charging time	≤ 5 hours
Working time	about 21 days (2 times for 2 mins)
Working temperature	5°C - 40°C
Waterproof class	IPX7

	The packaging of the device can be recycled according to local regulations.
	Symbol for corrugated cardboard (packaging material)

PARTS NAMES

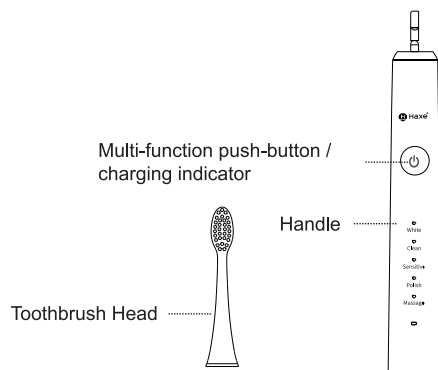TIP ATTACHMENT METHOD

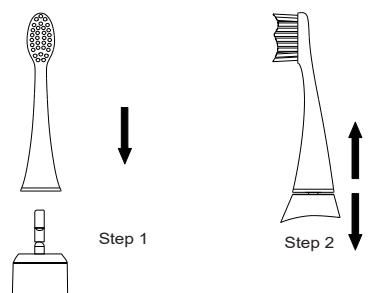FEATURES

1. The Haxe toothbrush ensures effective cleaning of the oral cavity, as it completes about 40 000 micro-brushing movements per minute.
2. The product offers 5 brushing modes, optimal for various users:
 - Whitening mode
 - Cleaning mode
 - Sensitive teeth mode
 - Polishing mode
 - Massage mode
3. The product is ergonomic, comfortable to hold in your hand. It does not load the wrist and does not slip away.
4. Its water-tightness rating IPX7 guarantees care-free use.
5. To charge the toothbrush, use the attached USB cable. Full loading is sufficient for about 21 days when brushing 2 times a day for 2 minutes.
6. After 30 seconds, change the brushing area. This will be signaled by a short pause. The device will automatically turn off after 2 minutes.
7. The memory function makes daily hygiene even simpler. The toothbrush will start in the last used mode.

Dear Customer!

Thank you for choosing Haxe products. We are committed to delivering sonic toothbrushes with an optimum quality/price ratio and superior performance to help our clients in their daily oral hygiene.

Read the user manual before first use. Keep the manual for future reference.

IMPORTANT SAFETY TIPS

1. Make sure that the mains voltage is compatible with the voltage on the device's nameplate.
2. Do not use the charging cord if it is damaged.
3. Keep the charging cord away from water. Do not immerse the cord in water or any other liquid.
4. Do not use the device if any of its parts is damaged. If the device needs to be repaired, follow the instructions in the "DEVICE FAILURE" section.
5. If the patient underwent jaw or gum surgery within two months prior to use, consult a dentist.
6. If severe bleeding occurs after using the toothbrush or if bleeding does not cease after one week of use, consult a dentist.
7. This toothbrush is a personal hygiene device and is not intended for use by many patients in a dental facility.
8. This device may be used by children over the age of 8 and by persons with limited physical and mental capacity only after they have received thorough instructions on safe use. Make sure that the user understands the risks and dangers associated with the use of the device.
9. Children must not clean and conserve the device unattended.
10. The product is not a toy, and children should not be playing with it.
11. The HX701 / HX702 toothbrush is designed to wash the tongue and gums. Do not use it for any other purposes.
12. Some people will feel discomfort (such as irritation or numbness) when using the toothbrush for the first time. This is normal when switching to a sonic toothbrush, and the process to get familiar with the toothbrush lasts about 1 week.
13. Protect the product from shocks and falls from height. If damage is noticed, immediately stop using the product.
14. The product can only be used with Haxe tip.
15. Do not use a tip with broken or crooked bristles, change the tip at least every 3 months.
16. If the user has a pacemaker or any other internal devices, they must consult their doctor in advance.
17. The product may only be used indoors. It should be protected from high temperatures and direct sunlight.

EXPLANATIONS OF THE SYMBOLS

	The device is marked with a crossed-out rubbish bin, in accordance with the European Directive 2012/19/WE on used electrical and electronic equipment (Waste Electrical and Electronic equipment - WEEE). Products marked with this symbol should not be disposed of or dumped with other household waste after a period of use. The user is obliged to get rid of used electrical and electronic equipment by delivering it to designated recycling point in which dangerous waste are being recycled.
	In accordance with directive 2013/56/UE on the disposal of batteries and accumulators, this product is marked with the symbol of a crossed-out waste bin. The symbol indicates that the batteries used in this product should not be disposed with regular household waste, but treated according to the legal directives and local regulations.

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